

ORTHO
SYSTEMS

2000
2026
Years in motion



Excellence in Neuro-Orthopedics

Your Patient – Our Contribution

Catalogue
Excerpt

Orthotics 14.0
Valid from July 2026

Stehfeld Standing Platform



Features

Sizes available for the **smallest to large** patients

Unrivalled **20-month guarantee**


Engraved **serial number** for traceability

Clearly illustrated instructions for use

+ Tilting connecting plate for **posture control exercises** during therapy

+ **Silicone bumper** to avoid damages and protect against injury

+ **Wheel locks** keep podium stationary

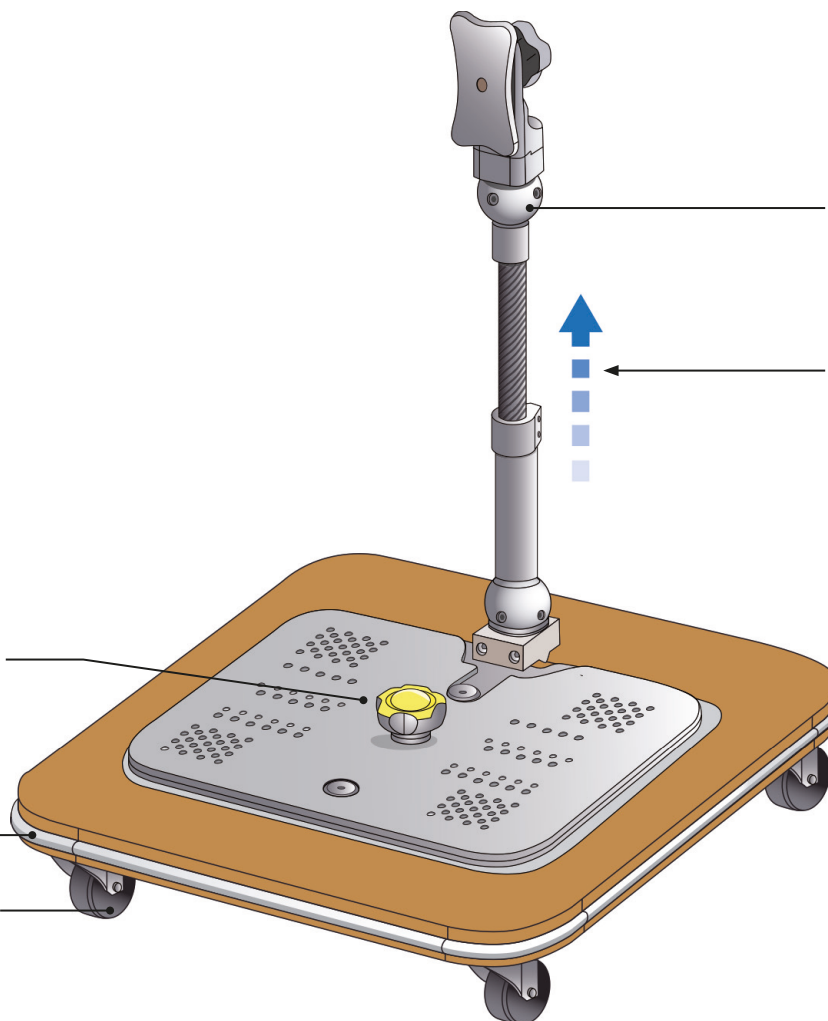

Standard

Fully assembled as one item for clear cost estimates and uncomplicated ordering


+ **Pyramid adapters** for precise adjustment, up to 10°


+ **Secure 3-point attachment** of the body shell at the feet and pelvic area

+ **Height adjustable** to accommodate growth

+ Unscrew **2 grips** to detach body shell and posture plate from the base


+ Applies the Helga Hirschfeld therapeutic principle to **induce postural reactions** by displacing the body and prompting a return to equilibrium

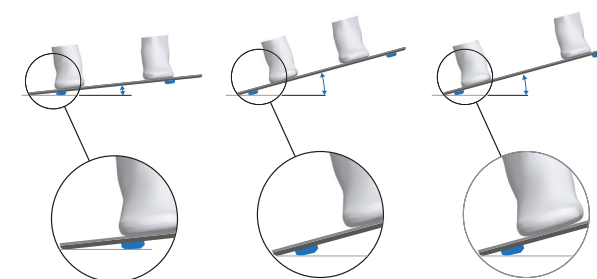

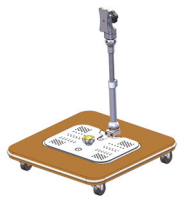


Article No.	Height Classification	Base (l x w)	Posture Plate	Foot Length	Foot to Foot (lateral)
Stehfeld-klein	65 cm to 100 cm	40 cm x 40 cm	22,5 cm x 29 cm	12,5 cm to 17 cm	17,0 cm to 23 cm
Stehfeld-mittel	80 cm to 130 cm	45 cm x 45 cm	25,5 cm x 36 cm	14,5 cm to 22 cm	20,5 cm to 31 cm
Stehfeld-groß	120 cm to 160 cm	50 cm x 50 cm	28,5 cm x 42 cm	16,5 cm to 25 cm	22,5 cm to 35 cm
Stehfeld-groß-breit	120 cm to 160 cm	50 cm x 70 cm	28,5 cm x 42 cm	16,5 cm to 25 cm	22,5 cm to 35 cm

With anti-tip protection to prevent toppling backwards

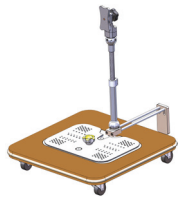


Article No. Combi-product	Height Classification	Base (l x w)	Posture Plate	Foot Length	Foot to Foot (lateral)
Stehfeld-KS-klein	65 cm to 100 cm	40 cm x 40 cm	22,5 cm x 29 cm	12,5 cm to 17 cm	17,0 cm to 23 cm
Stehfeld-KS-mittel	80 cm to 130 cm	45 cm x 45 cm	25,5 cm x 36 cm	14,5 cm to 22 cm	20,5 cm to 31 cm
Stehfeld-KS-groß	120 cm to 160 cm	50 cm x 50 cm	28,5 cm x 42 cm	16,5 cm to 25 cm	22,5 cm to 35 cm
Stehfeld-KS-groß-breit	120 cm to 160 cm	50 cm x 70 cm	28,5 cm x 42 cm	16,5 cm to 25 cm	22,5 cm to 35 cm

With anti-tip protection

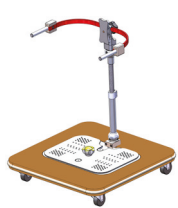

Article No. Combi-product	Height Classification	Base (l x w)	Posture Plate	Foot Length	Foot to Foot (lateral)
Stehfeld-HG-klein	65 cm to 100 cm	40 cm x 40 cm	22,5 cm x 29 cm	12,5 cm to 17 cm	17,0 cm to 23 cm
Stehfeld-HG-mittel	80 cm to 130 cm	45 cm x 45 cm	25,5 cm x 36 cm	14,5 cm to 22 cm	20,5 cm to 31 cm
Stehfeld-HG-groß	120 cm to 160 cm	50 cm x 50 cm	28,5 cm x 42 cm	16,5 cm to 25 cm	22,5 cm to 35 cm
Stehfeld-HG-groß-breit	120 cm to 160 cm	50 cm x 70 cm	28,5 cm x 42 cm	16,5 cm to 25 cm	22,5 cm to 35 cm

With hand grips and anti-tip protection

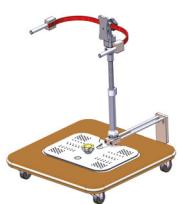


Article No. Combi-product	Height Classification	Base (l x w)	Posture Plate	Foot Length	Foot to Foot (lateral)
Stehfeld-HG-KS-klein	65 cm to 100 cm	40 cm x 40 cm	22,5 cm x 29 cm	12,5 cm to 17 cm	17,0 cm to 23 cm
Stehfeld-HG-KS-mittel	80 cm to 130 cm	45 cm x 45 cm	25,5 cm x 36 cm	14,5 cm to 22 cm	20,5 cm to 31 cm
Stehfeld-HG-KS-groß	120 cm to 160 cm	50 cm x 50 cm	28,5 cm x 42 cm	16,5 cm to 25 cm	22,5 cm to 35 cm
Stehfeld-HG-KS-groß-breit	120 cm to 160 cm	50 cm x 70 cm	28,5 cm x 42 cm	16,5 cm to 25 cm	22,5 cm to 35 cm

Body Shell with Posture Plate


Custom body shell used in conjunction with the Stehbrett supports primary deformities and selectively counters any possible secondary deformities.

Only a precise footbed and proper support of the ankle joint enable correct positioning of the knee and hip joints.


Amount of tilt can be individually adapted to the user by positioning supports on the undersurface.

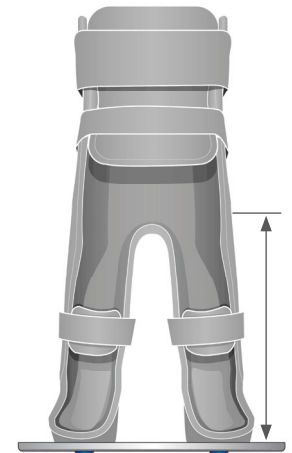


Minimum length from connecting plate to 4-point grip:
30 cm (small)
36 cm (medium and large)

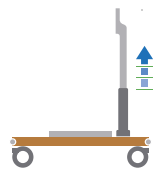


Growable

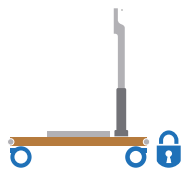


Stability

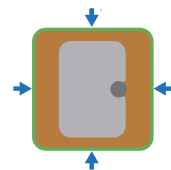


Protective Bumper

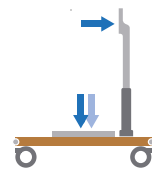


3-Point Fixation

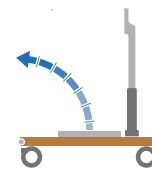


Removable

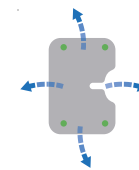

Postural Control Exercises

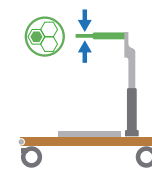

HG = Hand Grips

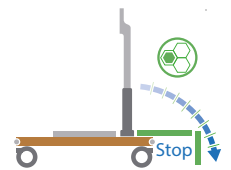

KS = Anti-tip Protection



Contact

You can reach us from
by telephone
per Fax
per E-Mail
Address
Contact persons

9 am to 5 pm (Germany)
+49 89 904 7575 0
+49 89 904 7575 5
kontakt@ortho-systems-mail.de
Merowingerstr. 14, 85551 Kirchheim, Germany
Jeff Matwichuk and Thomas Böckh

www.ortho-systems-inmotion.com