

ORTHO
SYSTEMS

2000
2026

Years in motion



Excellence in Neuro-Orthopedics

Your Patient – Our Contribution

Catalogue
Excerpt

Orthotics 14.0
Valid from July 2026

Modos 3D Hip Joint



Features

Sizes available for **small to large** patients

+ **Locking button** preloads the lock; engages automatically on extension.

+ **Unlocking button** pre-selects release; hip extension relieves pin tension and unlocks joint for sitting

+ **Pre-Selected** locking mechanism with 2-stage opening and closing


+ **Adduction is blocked.** 12° flexion stop can be increased as required for step length.

+ **Pelvic rotation** for efficient and energy saving gait. Joint axis permits **12° of internal and 12° of external rotation** independent of flexion/extension

+ **12° of abduction** increases stance stability and makes sitting more comfortable

Joint flexes **to 130°** when unlocked for sitting

Lower section made of high-strength, **hardened steel**

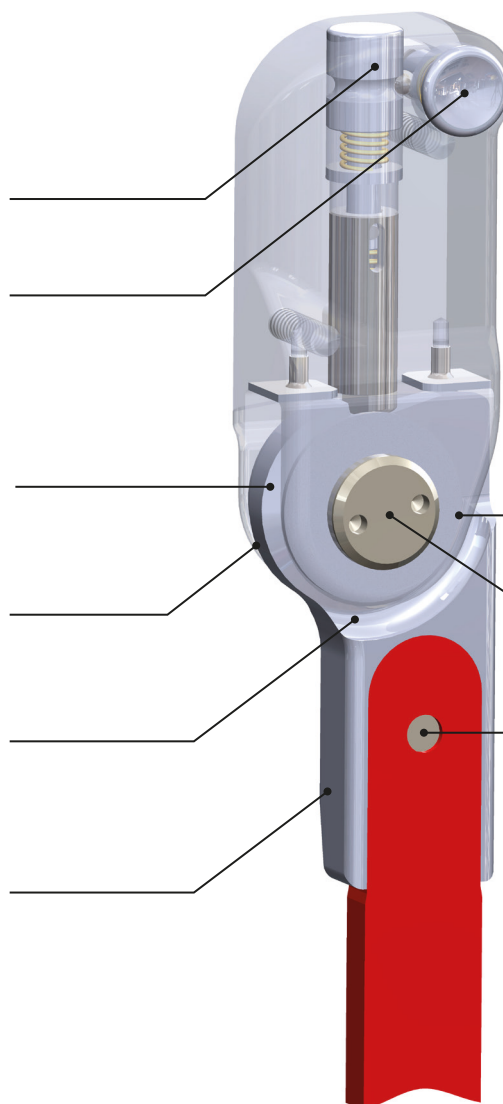

Standard

Fully assembled for clear cost estimates and uncomplicated ordering

Secure and adjustable connection to pelvic band

Upper section fabricated from **high-strength aluminium**

High-strength steel washers for smooth operation free of play

Maintenance-free, high-strength steel joint assembly with self-lubricating bearing


Machined upright and offset connection ensure **optimal transfer of force**

Unrivalled **20-month guarantee**

Engraved **serial number** for traceability

Clearly illustrated instructions for use

Fig. left

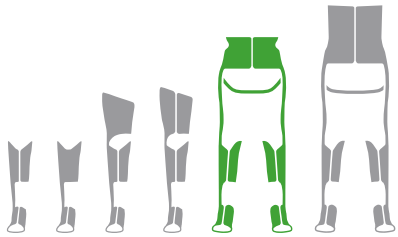


We recommend using our configurator to make your selection.
The configurator considers, among other factors, weight to height ratios.



Weight Classification (body weight)			Joint		
Article No.	Article No.	HKAFO	Width	Thickness	Height
Modos-3D-klein	Modos-3D-AG-klein	up to 20 kg	32,5 mm	23,0 mm	112,5 mm
Modos-3D-klein-Plus	Modos-3D-AG-klein-Plus	up to 30 kg	32,5 mm	23,0 mm	114,0 mm
Modos-3D-mittel	Modos-3D-AG-mittel	not yet available	-	-	-
Modos-3D-groß	Modos-3D-AG-groß	not yet available	-	-	-

AG = custom fabricated

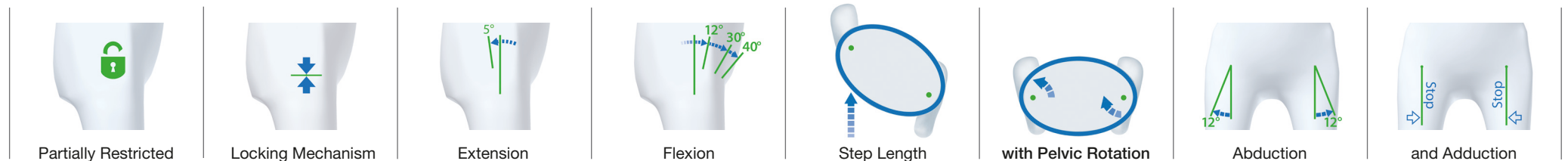



Contact

You can reach us from
by telephone
per Fax
per E-Mail
Address
Contact persons

9 am to 5 pm (Germany)
+49 89 904 7575 0
+49 89 904 7575 5
kontakt@ortho-systems-mail.de
Merowingerstr. 14, 85551 Kirchheim, Germany
Jeff Matwichuk and Thomas Böckh

www.ortho-systems-inmotion.com