

ORTHO
SYSTEMS

2000
2026
Years in motion



Excellence in Neuro-Orthopedics

Your Patient – Our Contribution

Catalogue
Excerpt

Orthotics 14.0
Valid from July 2026

Indi Knee Joint



Features

Waterproof for all situations

Sizes available for the **smallest to largest** patients

Available with **long or short uprights**

+ Maximum **extension 5°** and **flexion 115°**

+ **Posterior offset** proportionally adjusted to joint size

+ **Noise dampener** with self-locking setscrew for higher patient acceptance

Additional functions with the **Strecker add-on**

Fabricated from **high-strength titanium**

Standard

Fully assembled for clear cost estimates and uncomplicated ordering

Machined uprights and offset connection ensure **optimal transfer of force**

Upper section **CNC-machined from one piece** of material for increased stability

Precise assembly calibrated to 1/100mm with self-lubricating washers for smooth operation free of play

Maintenance-free stainless steel joint assembly with bronze bushing

Unrivalled **20-month guarantee**

Engraved **serial number** for traceability

Clearly illustrated instructions for use



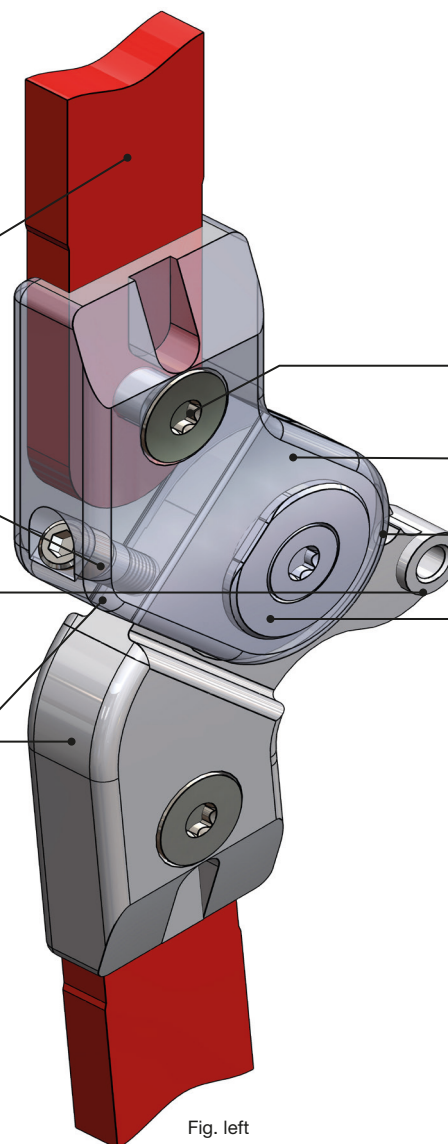



Fig. left

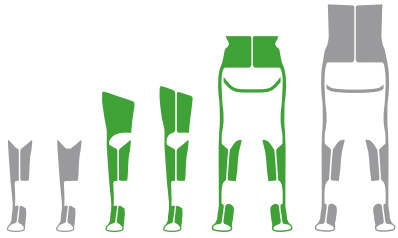

We recommend using our configurator to make your selection.
The configurator considers, among other factors, weight to height ratios.

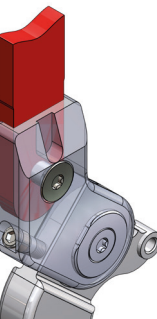

Weight Classification (body weight)				Joint				
Article No.	KAFO	HKAFO		Length	Width	Thickness	Weight	Posterior Offset
Indi-25-Plus	up to 30 kg	up to 25 kg		93 mm	25 mm	12 mm	83 g	15,5 mm
Indi-28	up to 40 kg	up to 35 kg		103 mm	28 mm	14 mm	119 g	16,5 mm
Indi-30	up to 50 kg	up to 45 kg		112 mm	30 mm	16 mm	152 g	17,5 mm
Indi-30-Plus	up to 70 kg	up to 60 kg		112 mm	30 mm	16 mm	152 g	17,5 mm
Indi-34	up to 85 kg	up to 85 kg		123 mm	34 mm	19 mm	230 g	22,0 mm

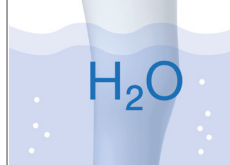								
Locked	Extension	Flexion	Posterior Offset	Modular	Force Options	No Grinding Necessary	Noise Dampening	Water Resistant



Contact

You can reach us from
by telephone
per Fax
per E-Mail
Address
Contact persons

9 am to 5 pm (Germany)
+49 89 904 7575 0
+49 89 904 7575 5
kontakt@ortho-systems-mail.de
Merowingerstr. 14, 85551 Kirchheim, Germany
Jeff Matwichuk and Thomas Böckh

www.ortho-systems-inmotion.com