

ORTHO
SYSTEMS

2000
2026
Years in motion



Excellence in Neuro-Orthopedics

Your Patient – Our Contribution

Catalogue
Excerpt

Orthotics 14.0
Valid from July 2026

Konso-0 Knee Joint



Features

Sizes available for **small to large** patients

Additional functions with the **Move-N** or **Move-S** add-ons

Joint flexes to **130°**

Available with a **long or short upright**

Ball plunger included
for temporary unlocking of joint

Fabricated from **high-strength aluminum**

+ Adjustable flexion angle from **0° to 40°** in 5° intervals

Available with a **contoured or straight lower bar**

Standard

Fully assembled for clear cost estimates
and uncomplicated ordering

Machined uprights and offset connection
ensure **optimal transfer of force**

Upper section **CNC-machined from one piece**
of material for increased stability

Precise assembly calibrated to 1/100mm with self-
lubricating washers for smooth operation free of play

Maintenance-free stainless steel joint assembly
with bronze bushing

Unrivalled **20-month guarantee**

Engraved **serial number** for traceability

Clearly illustrated instructions for use




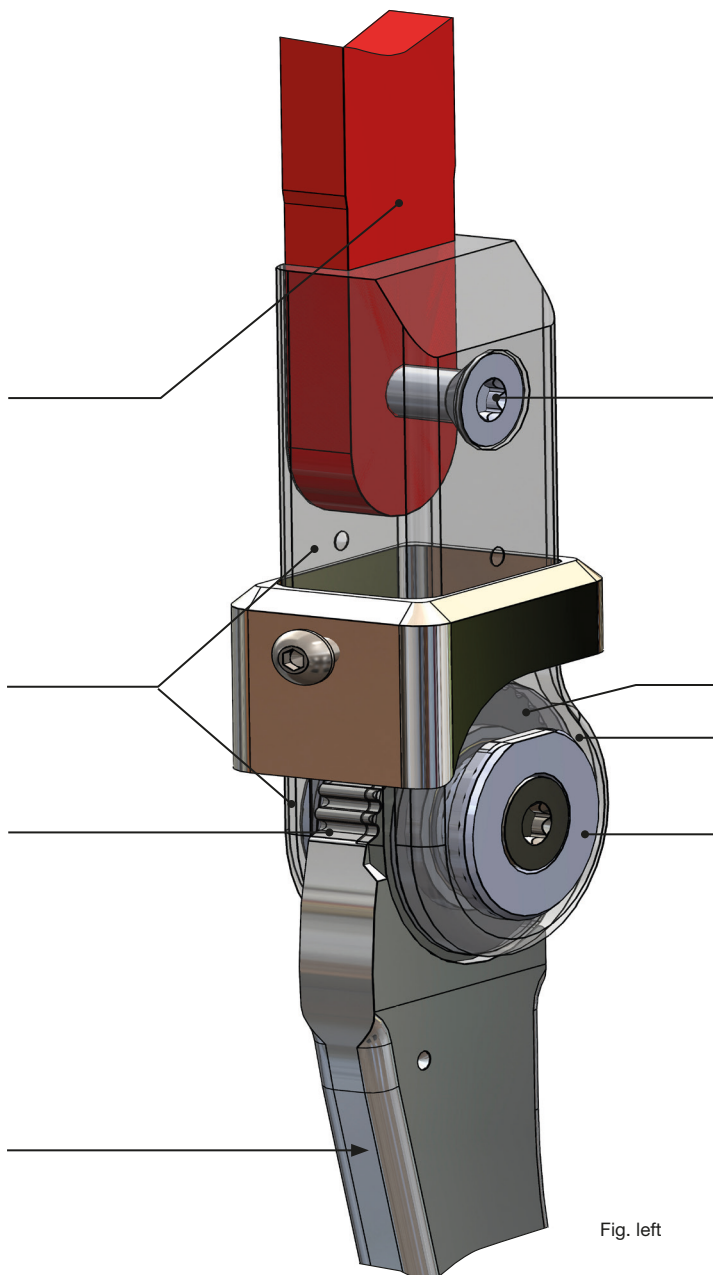


Fig. left

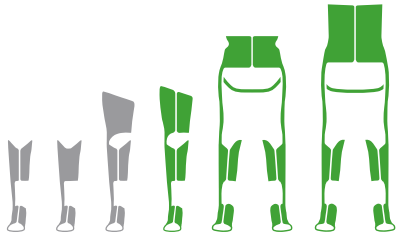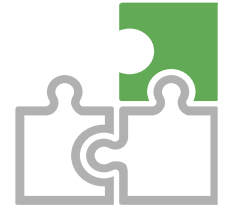



We recommend using our configurator to make your selection.
The configurator considers, among other factors, weight to height ratios.

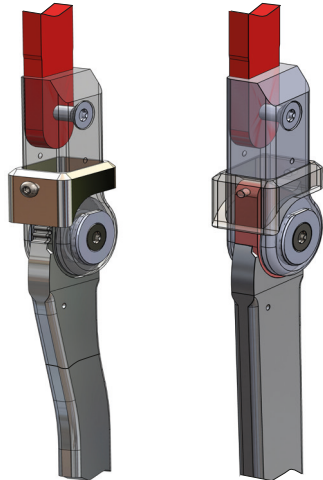

Weight Classification (body weight)			Joint			
Article No.	KAFO	HKAFO/RGO		Width	Thickness	Thickness at Drop Lock
Konso-0-25	up to 20 kg	up to 15 kg		25 mm	14,5 mm	17,5 mm
Konso-0-28	up to 30 kg	up to 20 kg		28 mm	16,5 mm	20,5 mm
Konso-0-30	up to 35 kg	up to 30 kg		30 mm	17,5 mm	21,5 mm
Konso-0-32	up to 50 kg	up to 40 kg		32 mm	19,5 mm	23,0 mm
Konso-0-36	up to 70 kg	up to 60 kg		36 mm	20,5 mm	24,0 mm
Konso-0-40	up to 85 kg	up to 85 kg		40 mm	22,0 mm	26,0 mm

Konso-0 contoured | Konso-0 straight

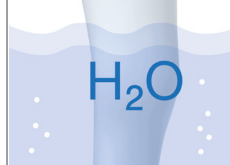								
Locked	Extension	Flexion	Posterior Offset	Modular	Force Options	No Grinding Necessary	Locking Mechanism	Water Resistant

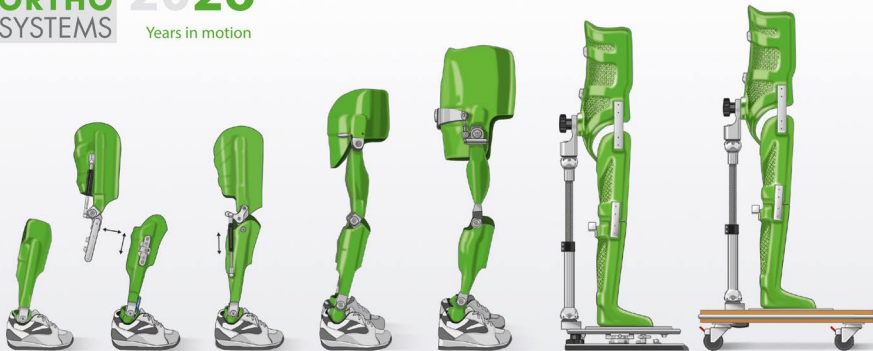

Contact

You can reach us from
by telephone
per Fax
per E-Mail
Address
Contact persons

9 am to 5 pm (Germany)
+49 89 904 7575 0
+49 89 904 7575 5
kontakt@ortho-systems-mail.de
Merowingerstr. 14, 85551 Kirchheim, Germany
Jeff Matwichuk and Thomas Böckh

www.ortho-systems-inmotion.com