

ORTHO
SYSTEMS

2000
2026
Years in motion



Excellence in Neuro-Orthopedics

Your Patient – Our Contribution

Catalogue
Excerpt

Orthotics 14.0
Valid from July 2026

Konso-5-Move-5-N



Features

- Integrated ball plunger** enables temporary unlocking of the joint
- + **Force unit** assists knee extension and helps when sitting down and standing up
- + **2-piece adjustable body** to vary flexion stop between **95°** and **115°**
- + **7 positions** to shift lever arm and increase or decrease torque
- Cushioned** flexion stop
- Wide range of **force options** available
- 2 different styles** of base plate for versatility: retrofit or incorporated

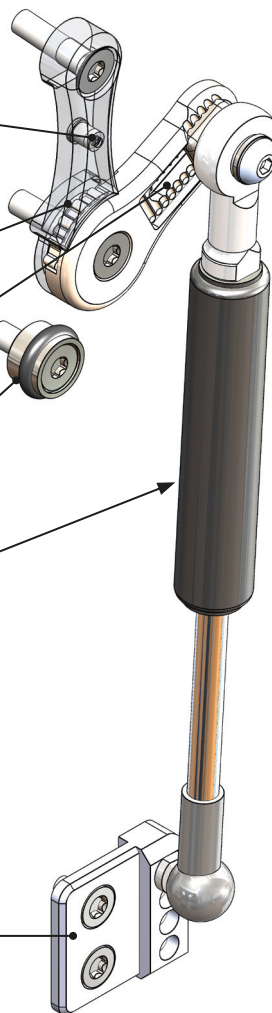


Fig. left


Standard

Fully assembled complete with everything needed for a functional add-on to the Konso knee joint


Available **separately** or as **combi-product** together with a Konso knee joint

Serial number for traceability

Clearly illustrated instructions for use

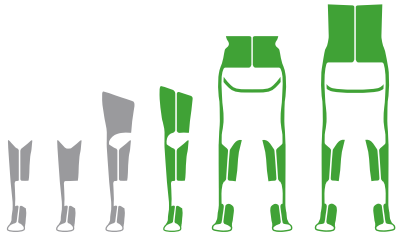

We recommend using our configurator to make your selection.
The configurator considers, among other factors, weight to height ratios.

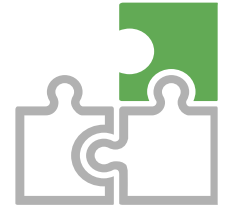


Weight Classification (body weight)

Article No. Combi-product	Weight Classification (body weight)		Article No.
	KAFO	HKAFO / RGO	
Konso-5-Move-5-N-25	up to 20 kg	up to 15 kg	Move-5-N-25
Konso-5-Move-5-N-28	up to 30 kg	up to 20 kg	Move-5-N-28
Konso-5-Move-5-N-30	up to 35 kg	up to 30 kg	Move-5-N-30
Konso-5-Move-5-N-32	up to 50 kg	up to 40 kg	Move-5-N-32
Konso-5-Move-5-N-36	up to 70 kg	up to 60 kg	Move-5-N-36
Konso-5-Move-5-N-40	up to 85 kg	up to 85 kg	Move-5-N-40

Konso-5 contoured | Konso-5 straight

Locked	Extension	Flexion	Posterior Offset	Modular	Force Options	No Grinding Necessary	Locking Mechanism	Water Resistant

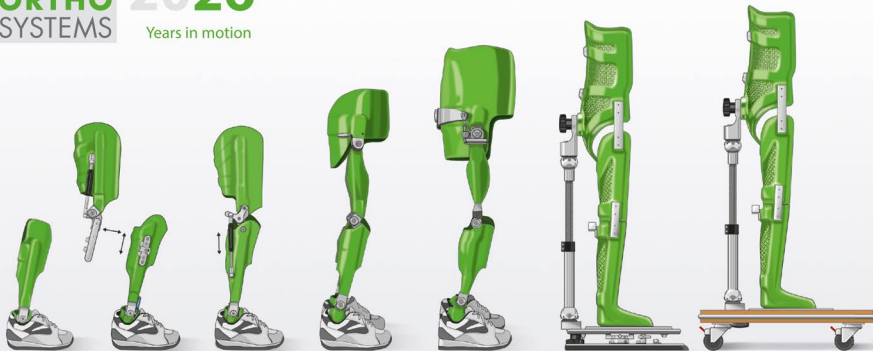

Contact

You can reach us from
by telephone
per Fax
per E-Mail
Address
Contact persons

9 am to 5 pm (Germany)
+49 89 904 7575 0
+49 89 904 7575 5
kontakt@ortho-systems-mail.de
Merowingerstr. 14, 85551 Kirchheim, Germany
Jeff Matwichuk and Thomas Böckh

www.ortho-systems-inmotion.com